

HL7 Output

Overview:

This guide provides instruction on configuring the Medical-Objects Download Client (Capricorn) to output HL7 Results to a set folder.

If you require any assistance with the installation or configuration, or if you have any questions, please contact the Medical-Objects Helpdesk on (07) 5456 6000.

If you have not installed the Medical-Objects Download client, Please see the guide for [Installing Capricorn on Windows](#).

- 1 [Overview:](#)
- 2 [Configuring Capricorn Software](#)
 - 2.1 [Accessing the Capricorn Configuration Window](#)
 - 2.2 [Capricorn Configuration Settings](#)

Configuring Capricorn Software

Accessing the Capricorn Configuration Window

Note:

In some cases, the Capricorn may be configured as a service. See [this guide](#) for instructions on how to launch the Capricorn if it is set up as a service.

1. To launch the HL7 Tray application, Go to **Start Menu, All Programs, Medical-Objects** and select the **Medical-Objects Capricorn**. Alternatively, You can go to **C:\MO\Capricorn** and run the **Capricorn.exe**.
2. The following icon will appear in the tray. Double click on it to bring up the log window.




3. When the log window is displayed, Click on the **Configuration** icon highlighted below, or select **Utility** and select **Configuration**.



Capricorn Configuration Settings

1. Click on **Integration Options**.

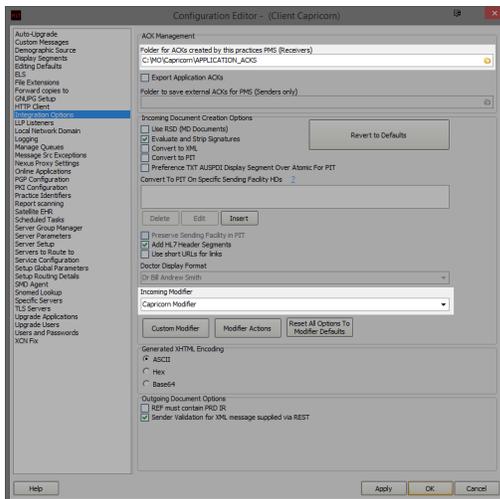


2. Set the **ACK Management** path to **C:\MO\Capricorn\APPLICATION_ACKS**.
3. Set the **Incoming Modifier** to **Capricorn Modifier** as shown above.



Alternatively the field can be left blank for a messages to be delivered as close to they were originally sent as possible.

4. Click on **Server Parameters** in the left pane.

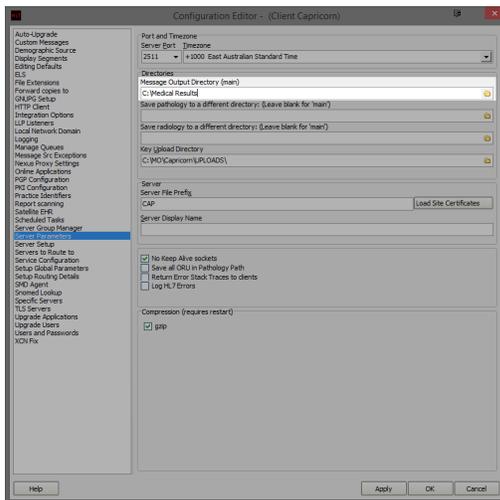

5. Under the **Message Output Directory (main)**, set the results import path.

The path will be a folder you want to use to store the HL7 results or messages that you receive.

For this installation we will just create a folder in the C: drive called **Medical Results**.
The file path in this case will be **C:\Medical Results**.

6. Press the **Apply** button when configurations are complete. All future messages will be delivered to this folder in HL7 format.