

PROMs User Guide

Overview

This article will cover how to create and use the Medical-Objects PROMs feature in a waiting room scenario. If you have any questions or require any further information, please contact the Helpdesk on (07) 5456 6000.

Requirements


PROMs is a premium add-on feature that is accessed directly from the Explorer Online web-based portal. This is designed to run on desktops, laptops and modern tablets running Windows, MacOS or Linux. Mobile phones are not currently supported.

Creating a Questionnaire

If you require additional patient questionnaires to be digitised or added, please contact the Sales team [here](#)

To create a PROMs questionnaire.

1. Click the **Documents Patient Questionnaire** (under Create) in the navigation bar at the top.

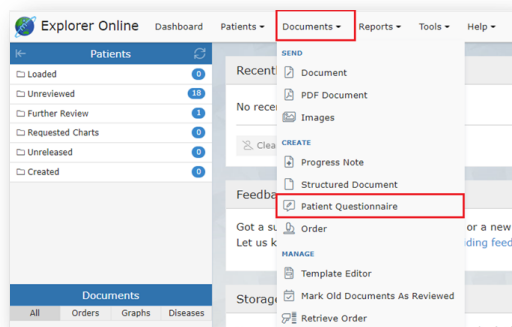

2. Fill in the form details. Add the patient's **family name**, **given name** and **date of birth**.
3. Choose the form template that you wish to use.
4. Select which provider the completed form will be sent to and press the **Add** button if you would like to send copies to additional providers.

Create Patient Questionnaire

FAMILY NAME:
TEST

GIVEN NAME: Patient DATE OF BIRTH: 01/06/2018

FORM TEMPLATE:
Calculate DASS 21 score

ADDRESSED TO:

+ Add X Clear


COPIES TO

Cancel Submit

- 1 Overview
- 2 Requirements
- 3 Creating a Questionnaire
- 4 Accessing a Questionnaire
 - 4.1 Using the QR Barcode
 - 4.2 Using the Link
- 5 Completing the Questionnaire
- 6 Viewing Completed Questionnaire
- 7 Cumulative Results
- 8 History of Created PROMs
 - 8.1 Understanding the Report

Accessing a Questionnaire

Once you click Submit, a **QR Code** and **link** option will appear.

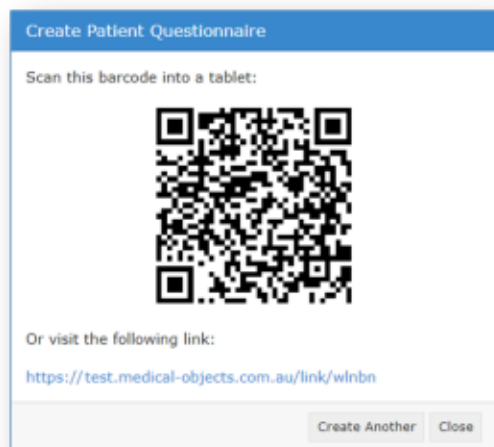

Using the QR Barcode

You can use the **barcode** option to transfer the form to another device, such as a tablet.

Open the device's **Camera app** or **QR code reader** and hold it over the **barcode**. Then, tap the notification to open the questionnaire. Next, hand the form to the patient to complete.

Using the Link


Clicking the Link will open the form in a new tab. You may also manually copy the link to be able to send it using alternate methods (e.g. via Email).

Please note that clicking on the link will bring you directly to the Patient Questionnaire. Only select the **link** option if you or another staff member wish to complete the form for your patient.

Completing the Questionnaire

The patient can answer the questions by typing into the required fields/selecting the required check boxes. There is also an option to:

- Reset the form
- Save it as a draft
- Submit it

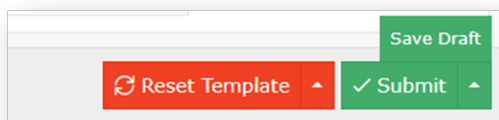


Once completed and submitted, the form automatically determines where the patient lies on the scale and an instant PROMs report is created. This report is then sent directly to Explorer Online or delivered to your compatible practice software.

Viewing Completed Questionnaire

To view the delivered PROMs report in Explorer Online, load the patient in the **sidebar**.

If the patient is not loaded in the **sidebar**, perform a patient lookup. (For more information on how to lookup patients, please click [here](#) to view the Explorer Online User Guide)

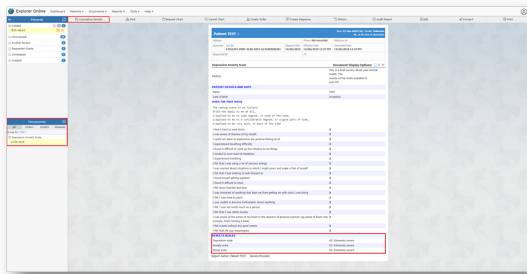

Once the patient is loaded, select their questionnaire report from the **Document** section of the sidebar.

Here, you can view the report details, the patient's responses, and the result of the total score.


Alternatively, if you are using compatible practice software the completed PROMs report can also be viewed there.

Cumulative Results

The **Cumulative Results** tab is another feature of Explorer Online that you can access when viewing report details. This feature allows you to compare the patient's total scores (numerical data only) after they have submitted the same survey multiple times.

If you would like to preview a document in the results area:


- Click on the date the survey was completed.
- Tick the **Total** checkbox or the specific **Test Name** to graph the changes in the patient's condition over time.

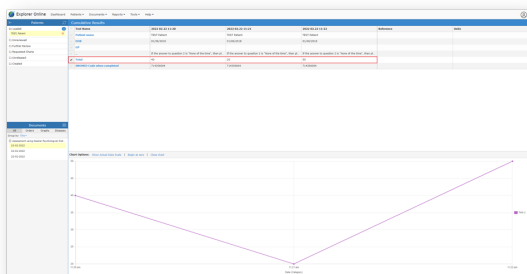


History of Created PROMs

To view the history of created PROMs, see if they have been completed or access the **QR Code** and **link** again;

1. Select **Patient Questionnaires** from the **Reports**

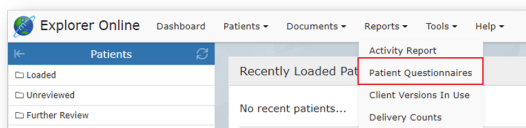

2. Choose from the following filter options

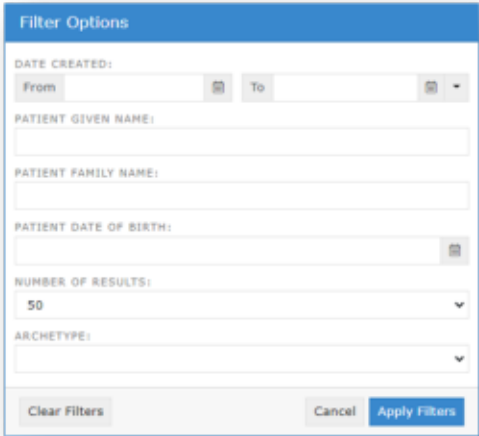

The 'Filter Options' dialog box contains the following fields and controls:

- DATE CREATED:** A date range selector with 'From' and 'To' input fields and calendar icons.
- PATIENT GIVEN NAME:** A text input field.
- PATIENT FAMILY NAME:** A text input field.
- PATIENT DATE OF BIRTH:** A date input field with a calendar icon.
- NUMBER OF RESULTS:** A dropdown menu currently set to '50'.
- ARCHETYPE:** A dropdown menu.
- Buttons at the bottom: 'Clear Filters', 'Cancel', and 'Apply Filters'.


- a. The number of results and the patient details.
 - b. Choose from the date created range
 - c. Select the form type
3. Once you are finished, hit the **Apply Filters** button.

Click the **plus button** to create a new PROMs questionnaire.

Select the **refresh button** to update results and click the **filters button** to edit the filter options.



You can also sort, hide, show or reset columns via the **grid menu**. Additionally, forms can be sorted into columns and ordered by **Patient**, **Date of Birth**, **Date Created**, **Archetype**, and **Date Satisfied**.

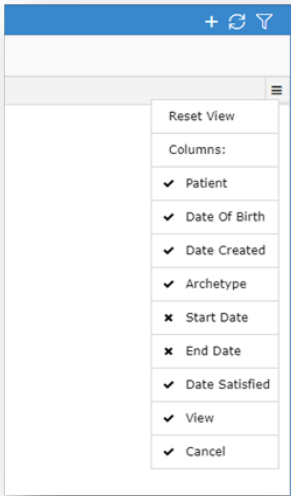

Understanding the Report

Patient	Patient name. SURNAME, First name
Date of Birth	Patient's Date of Birth
Date Created	The date the document was created
Archetype	The title of the Patient Questionnaire that was created

Date Satisfied	The date the Questionnaire was completed
View	View the QR code and link of Questionnaire
Cancel	Can cancel the request to complete the Questionnaire
<div>Date Created ▼</div>	Click the top of each column to sort by ascending or descending values.