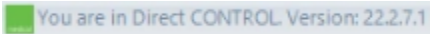


Direct CONTROL

Overview

This article is to assist with the setup of Medical Objects and Direct Control

Requirements

1. Direct CONTROL version 22.2.7.1 or higher is required. This can be checked in the top right corner of Direct CONTROL
You are in Direct CONTROL Version: 22.2.7.1
2. Capricorn 3.2.13.138 or higher is required.

Configuring Capricorn Software

Accessing the Capricorn Configuration Window

In some cases, the Capricorn may be configured as a service. See [this guide](#) for instructions on how to launch the Capricorn if it is set up as a service.

1. To launch the HL7 Tray application, Go to **Start Menu, All Programs, Medical-Objects** and select the **Medical-Objects Capricorn**. Alternatively, You can go to **C:\MO\Capricorn** and run the **Capricorn.exe**.
2. The following icon will appear in the tray. Double click on it to bring up the log window.




3. When the log window is displayed, Click on the **Configuration** icon highlighted below, or select **Utility** and select **Configuration**.



Capricorn Configuration Settings


1. Click on **Ack Management**. Set the **Folder for Acks created by this practices PMS (Receiver)** to the HL7 Out path in the walkthrough here: [HL7Out](#)

- 1 [Overview](#)
- 2 [Requirements](#)
- 3 [Configuring Capricorn Software](#)
 - 3.1 [Accessing the Capricorn Configuration Window](#)
 - 3.2 [Capricorn Configuration Settings](#)
- 4 [Direct CONTROL Import Configuration](#)
 - 4.1 [HL7 In](#)
 - 4.2 [HL7 Out](#)
- 5 [Inbox](#)

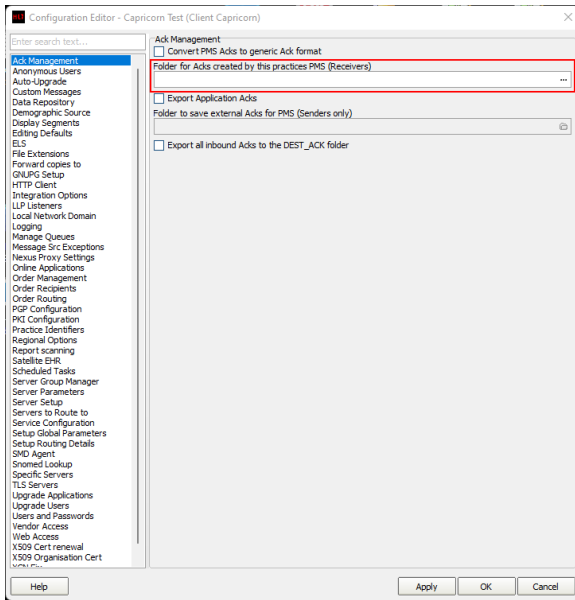


2. Click on **Integration Options**.

3.1 Check that the **Incoming Modifier** is set to **Direct CONTROL**

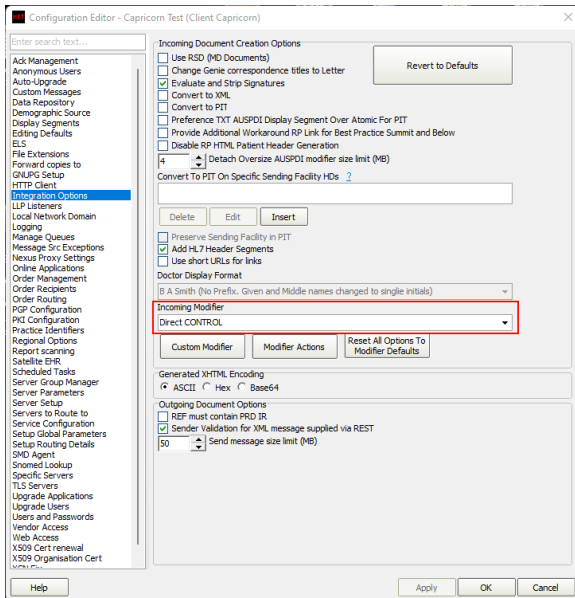

3. Select **Server Parameters** from the left panel.

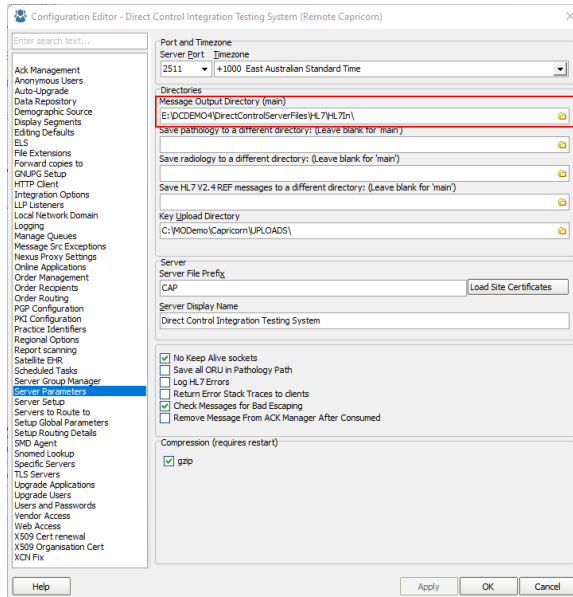

The paths above are set as follows:

```
(Capricorn installed on a Direct CONTROL Server Machine) --> E:\
\InstallFolder\DirectControlServerFiles\HL7\HL7In\
```

```
(Capricorn installed on a Direct CONTROL Client Machine) -->
\\server\InstallFolder\DirectControlServerFiles\HL7\HL7In\
```


If the folder does not exist, You must create it.

4. The Medical-Objects configuration is complete. Click **Apply** and **OK**.

Direct CONTROL Import Configuration

In Direct CONTROL navigate to Utilities Configure

HL7 In

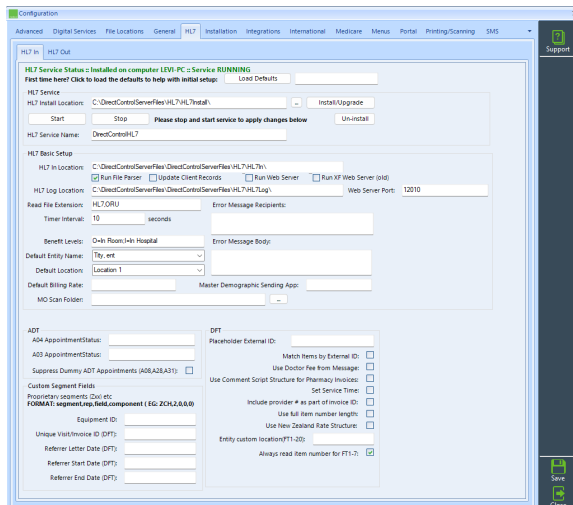

If paths are not populated click "Load Defaults" then "Install/Upgrade"


Change settings as required.

Set the download path in **HL7 In Location** as the first setting in **HL7 Basic Setup**.

Select "Start" button

HL7 Out

This is currently only for return acknowledgements sending messages directly from Direct CONTROL has not been finalised.

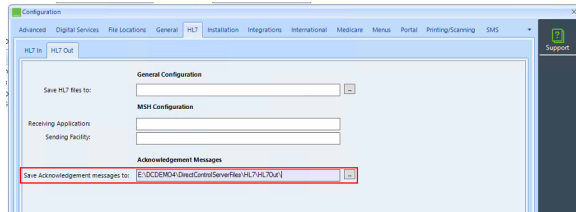


Set the "Save Acknowledgement messages to:" path.

Inbox

Regular Tasks --> Communications/Results

If Direct CONTROL was not able to find a match for an existing Client, the documents will appear in the **Unallocated** tab.

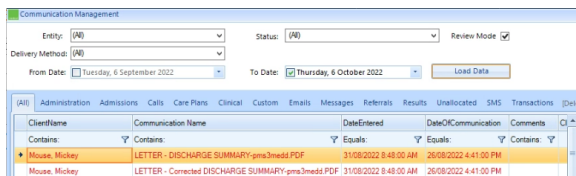

ClientName	Communication Name	DateEntered	DateOfCommunication	Comments
Mouse, Mickey	LETTER - DISCHARGE SUMMARY-pms3medd PDF	31/08/2022 8:48:00 AM	26/08/2022 4:41:00 PM	
Mouse, Mickey	LETTER - Corrected DISCHARGE SUMMARY-pms3medd PDF	31/08/2022 8:48:00 AM	26/08/2022 4:41:00 PM	