

Editing and sending Typist Messages

Overview


- 1 Overview
- 2 Editing and Sending Typist Messages

This guide provides a brief overview on how to Edit and send messages saved by Typists from within the Medical-Objects Explorer (Meridian).


If you require any assistance with the installation or configuration, or if you have any questions, please contact the Medical-Objects Helpdesk on (07) 5456 6000.

Editing and Sending Typist Messages


1. Open your Medical Objects software and login.
2. Select the **Note Pad** from the button toolbar at the top of the Explorer indicated in the screen shot below.

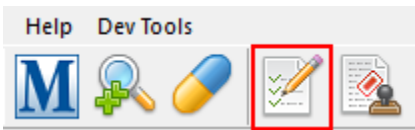


3. If there are any messages to edit they will show under the **"For Editing"** title in the left hand menu indicated in the screen shot below.

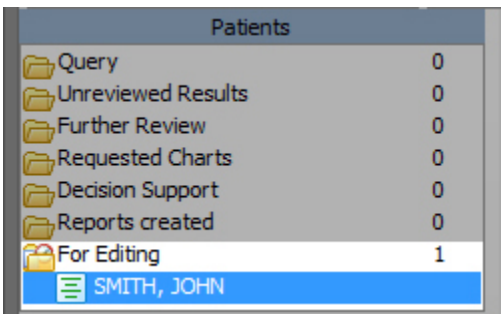

4. Next click on **Edit report** in the report to make any changes you require.

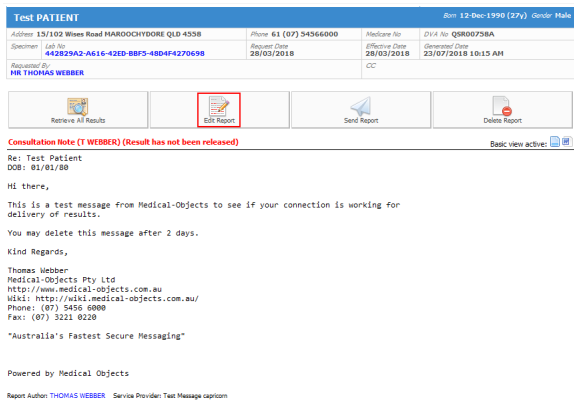

5. Make any changes to the message you require in the editor that then opens up. When done click on the **Mail** icon under the file menu indicated in the screen shot Below.

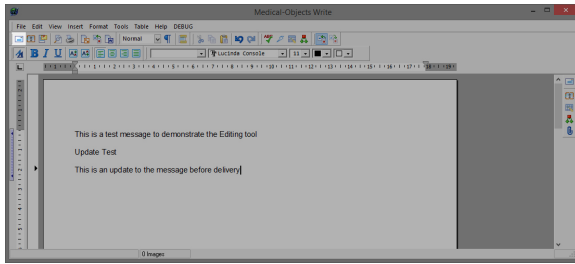

- In the following screen the minimum information required is Patient Name, who the report is from and an address for who it is going to. To address the report, click on the **Magnifying Glass** in the top right corner of the screen and use the consequent window to find your recipient. When you are ready to send the result make sure that **"For Approval before delivered"** is no longer ticked and you will notice that the **Save** button turns then turns into an **OK** button.

- You will then need to click **send** for the result to go to the recipients.

Note:

If there is nothing that you need to change in the report you can click **"Send Report"** without clicking **"Edit Report"** button when no editing is needed.

Test Patient

Address: 15/192 Miles Road HARBOURSIDE QLD 4038 Phone: 61 (07) 34364000 Medicare No: 010000738A DOB: 12-Dec-1990 (27y) Gender: Male

Specimen: Lab No: 442829A2-A616-42ED-8BF5-48DAF4270698 Request Date: 28/03/2018 Effective Date: 28/03/2018 Generated Date: 23/07/2018 10:15 AM

Reported by: MR THOMAS WEBBER

Buttons: Retrieve All Results, Edit Report, **Send Report**, Delete Report

Consultation Note (T WEBBER) (Result has not been released)

Re: Test Patient
DOB: 01/01/80

Hi there,

This is a test message from Medical-Objects to see if your connection is working for delivery of results.

You may delete this message after 2 days.

Kind Regards,

Thomas Webber
Medical-Objects Pty Ltd
http://www.medical-objects.com.au
Wkiki: http://wiki.medical-objects.com.au/
Phone: (07) 3436 4000
Fax: (07) 3221 8220

"Australia's Fastest Secure Messaging"

Powered by Medical Objects

Report Author: THOMAS WEBBER Service Provider: Test Message caption

- Finally you will need to check to see if the report has made it to the destination. Please see the guide [Checking Results Delivered](#) for further information. You can also see the [Viewing Results Marked for Further Review](#) section for more information about documents marked for Approval.