

# Monet

## Overview:

This guide provides instruction on configuring both Medical-Objects Capricorn and Monet software to import electronic results.

If you require any assistance with the installation or configuration, or if you have any questions, please contact the Medical-Objects Helpdesk on (07) 5456 6000.

If you have not Installed the Medical-Objects Download client, Please see the guide for [Installing Capricorn on Windows](#).

- 1 [Overview:](#)
- 2 [Configuring Capricorn Software](#)
  - 2.1 [Accessing the Capricorn Configuration Window](#)
  - 2.2 [Capricorn Configuration Settings](#)
- 3 [Configuring Monet](#)
  - 3.1 [Importing Pathology & Radiology Results](#)

## Configuring Capricorn Software

### Accessing the Capricorn Configuration Window

Note:



In some cases, the Capricorn may be configured as a service. See [this guide](#) for instructions on how to launch the Capricorn if it is set up as a service.

1. To launch the HL7 Tray application, Go to **Start Menu, All Programs, Medical-Objects** and select the **Medical-Objects Capricorn**. Alternatively, You can go to **C:\MO\Capricorn** and run the **Capricorn.exe**.
2. The following icon will appear in the tray. Double click on it to bring up the log window.




3. When the log window is displayed, Click on the **Configuration** icon highlighted below, or select **Utility** and select **Configuration**.




### Capricorn Configuration Settings

1. Click on **Ack Management**.
2. Set the Ack path to **C:\MO\Capricorn\APPLICATION\_ACKS**

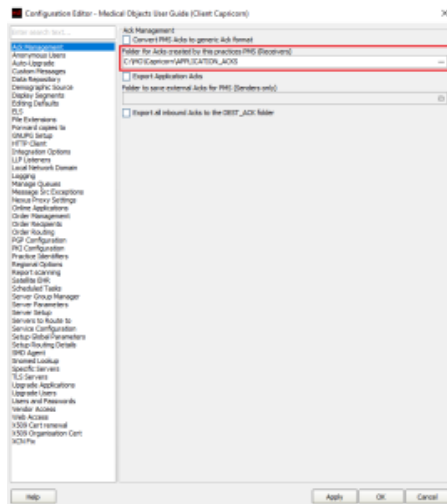


### 3. Click on **Integration Options**.

#### 3.1 Check that the **Incoming Modifier** is set to **Monet**

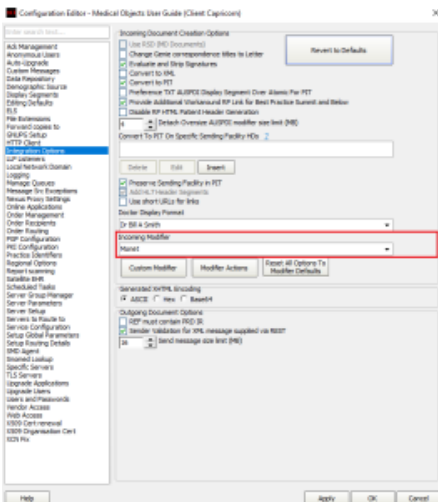

### 4. Select **Server Parameters** from the left hand panel.

### 5. Enter the path to drop off results to - for this exampled we'll use **C:\Pit\Hold**

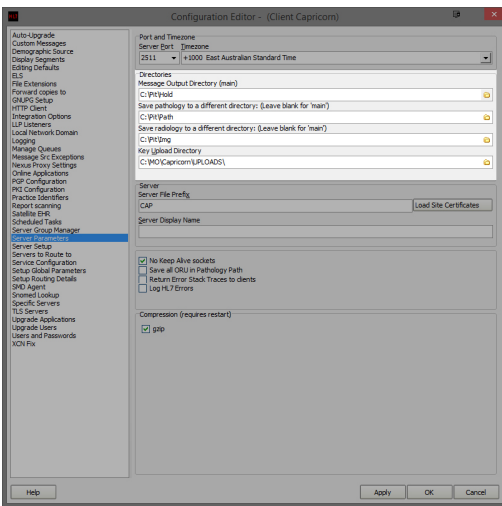

If receiving from Pathology and Radiology, Set the file paths as follows.

Pathology --> **C:\Pit\Path**

Radiology --> **C:\Pit\Img**


Note:

If the folders do not exist you must create them.


6. Click **Apply** and **OK** to save your configuration.

## Configuring Monet


1. To access the result input configurations go to the '**Monet palette**' window of Monet.
2. Click on the **Check Result** Icon in the left pane.

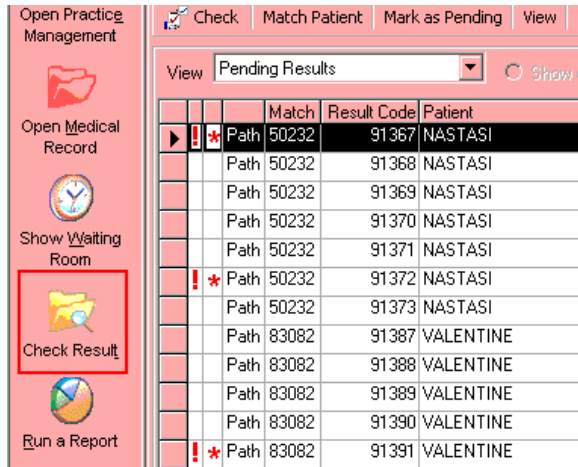

3. Click on the **Configure** Button shown below.

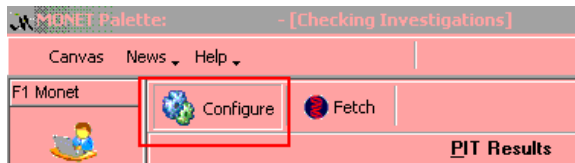


4. The results input path screen will be displayed.

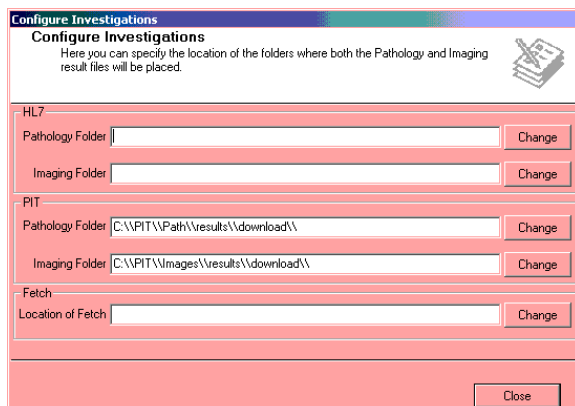


## Importing Pathology & Radiology Results


1. Open **Monet**, and from the **Palette** screen, click the **Check Results** button to open the **Check Results** screen.

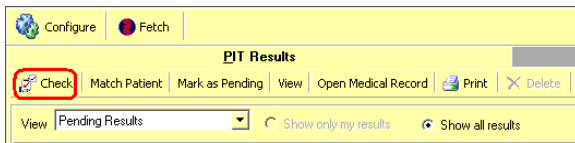

2. Click the **PIT** tab and then click the **Check** button to import the new results into Monet.
3. The **PIT Pathology Results** message box will appear advising that the results have been imported. Click the **OK** button.

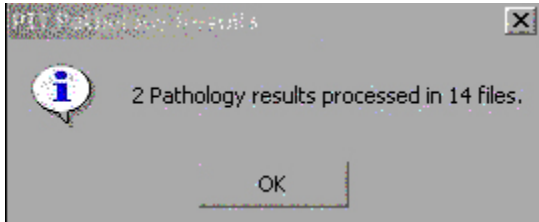

4. Repeat steps 1-3 for each Pathology and Radiology company from the external downloading software.
5. Each pathology & radiology company data must be downloaded and subsequently imported into the **Check Results** screen one company's group of results at a time. Do not combine the process of downloading the data for each company and then attempt to import them together.
6. The imported results will display in the **Check Results** screen. The results will automatically match to a Patient if all the Patient's details within the result file (the . pit file) match the Patient's demographics in Monet.
7. Where the imported results cannot match to a Patient's file, a red cross will be displayed in the **Match** column.

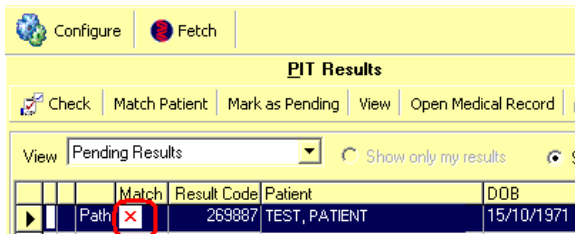

8. Click the **Match** button on the top toolbar to display the **Patient Match** screen.
9. Type in the Patient's name to perform the search and then click the **Match Again** button.

