

How to setup a basic send queue

Overview


This guide provides an overview on setting up a basic sending queue within the Medical-Objects Eclipse software.

If you have any questions or require any further information, Please contact the Medical-Objects Helpdesk on (07) 5456 6000.


- 1 [Overview](#)
- 2 [Setting up a send Queue](#)
- 3 [Running as a service](#)

Setting up a send Queue


1. Enter the configuration and go to **"Client Logging"** Set the Admin e-mail address where you would like notifications of problems to be delivered to.

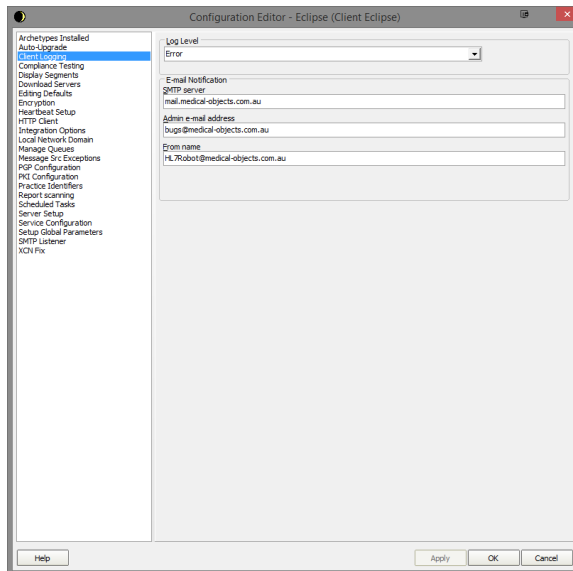


2. Go to **"Heartbeat Setup."**
3. Select the **Active** tick box.
4. Set interval to **30 minutes**.
5. Select the **"Server to Upload to"** from the drop down and choose one of the servers.

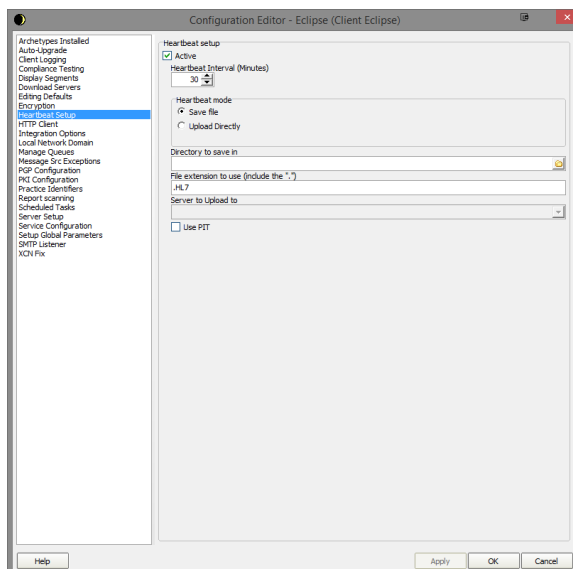

6. In **"Scheduled Tasks"** choose **"Zip Archive Message Store"**.
7. Choose the time that you would like to perform this task.
8. Set the **"File age to Archive (Days)"** to **10 days** if not less.

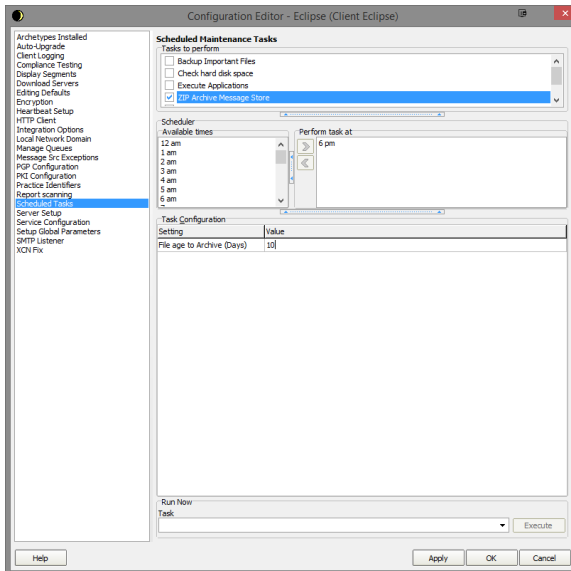

9. Choose **"Setup Global Parameters"**.
10. Make sure **"Seconds to wait before processing files"** is set to at least **3 seconds**.
11. This is to ensure that the file dumped is no longer being written before it is processed so if it can be set as high at **10 seconds** it should decrease issues.

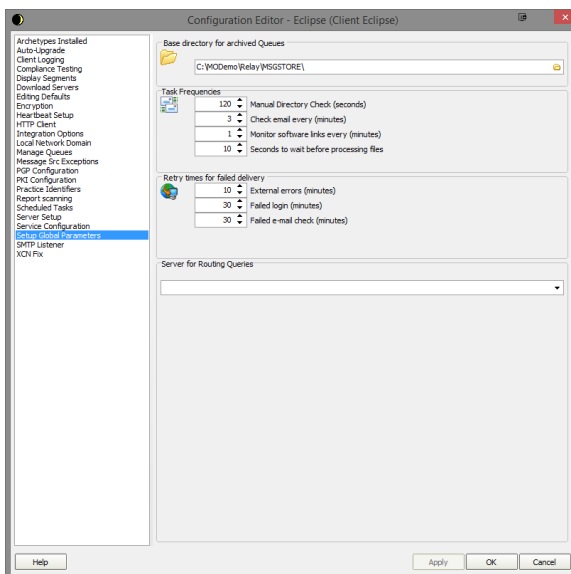


12. Choose **"Manage Queues"**.
13. Click the **"Add"** button.
14. Enter a name for the queue.
15. Enter the folder where the files are to be picked up from.

Note:




You can tick and untick the Active button at the bottom to disable and enable the queue.


16. Finally click the **"Setup"** button for the next step.

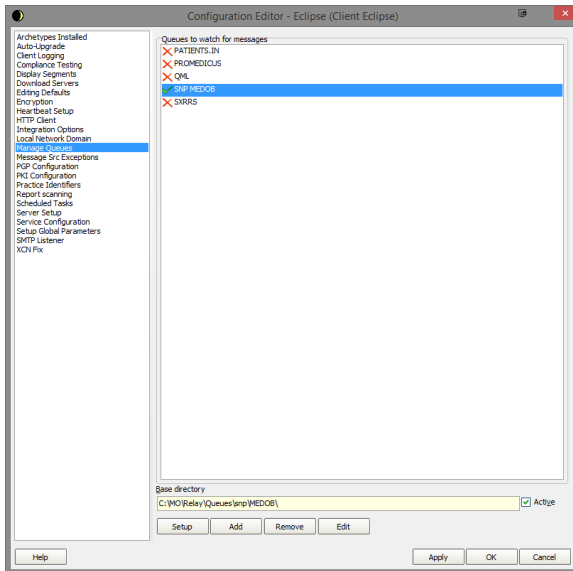


17. Click on the “Add” button.
18. Add each of the modifiers as per the screenshot in step 20.
19. For the “Store Message in Folder” choose a folder to keep the original file in for situations when we may need to trace any changes to the file from start to finish.
20. Make sure the tick boxes are selected as per the screenshot.

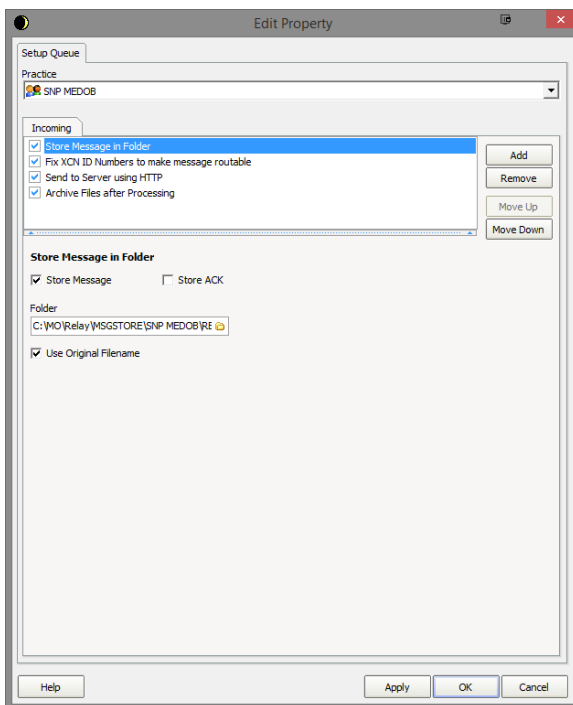

21. In the “Send to Server using HTTP” modifier set as below. Choose from the “Server name” one of the servers listed in the drop down.

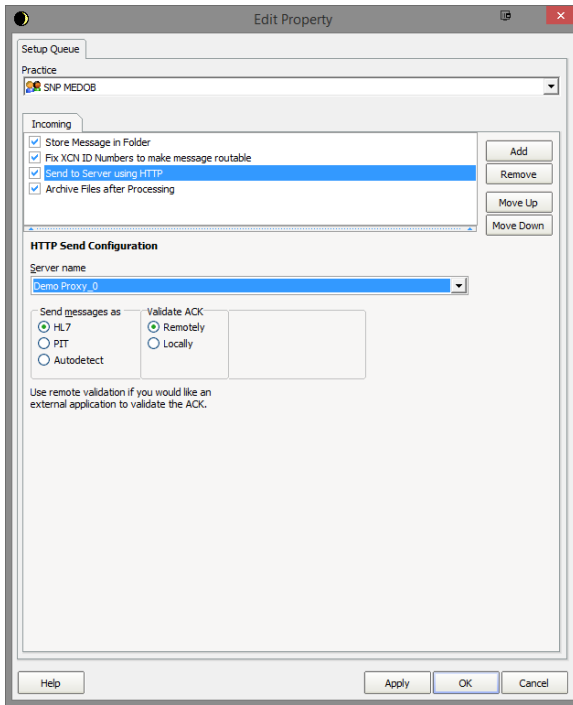


22. In “Archive Files after Processing” make selections as per the screenshot below.

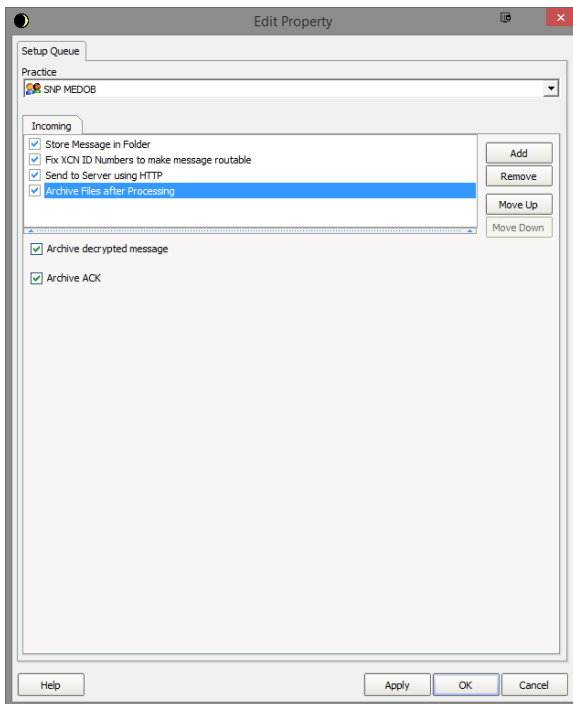

Running as a service

1. Choose “Service Configuration”.

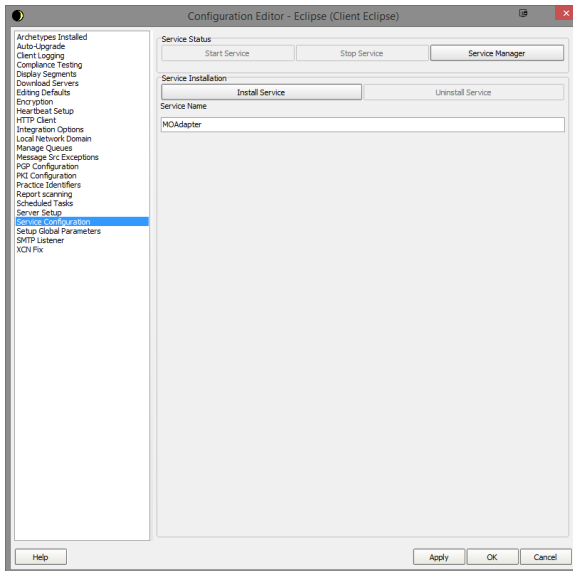

2. Before you press Install Service you can set the Service Name. Generally this is set as either **MOAdapter** or **Eclipse**.
3. When ready select Install Service and it should install into your services.
4. Exit the configuration and shutdown the eclipse and you can start up the service. **Do not run the tray and service at the same time or they will conflict trying to access the same configuration files.**

Note:



The Service by default installs as a local system user and you may have to change this user if it can't process files from the share folders it accesses.