

ZedMed Sending

Overview

This guide provides instruction on configuring both Medical-Objects Capricorn and ZedMed software for sending. This step is designed to be included in the Capricorn configuration for ZedMed [here](#).

If you require any assistance with the installation or configuration, or if you have any questions, please contact the Medical-Objects Helpdesk on (07) 5456 6000.

If you have not Installed the Medical-Objects Download client, Please see the guide for [Installing Capricorn on Windows](#).


- 1 [Overview](#)
- 2 [Sending the referral](#)
- 3 [Sending Batch Letters via ZedMed](#)

Sending the referral


1. Open up the patient you wish to send a report in regards to. Next click on the Start Encounter button within Zedmed or press F10 on your keyboard.



2. Click on the **Referral** button or press **CTRL+R** on your keyboard.

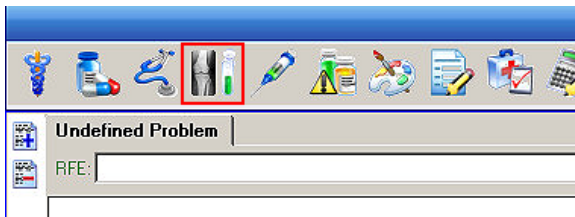

3. Select the **New Referral** tab, and select **Referral Letter**.

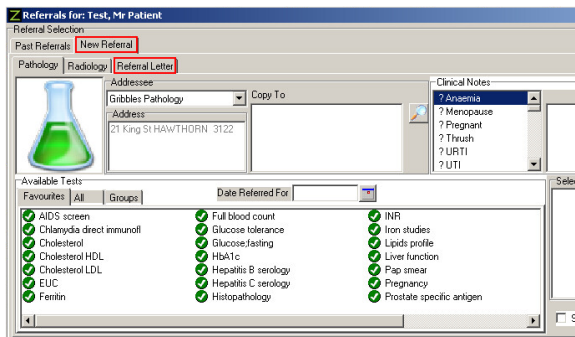

4. Under the **Addressee** field, Click the Magnifying Glass to search for the doctor you wish to send to.

Referrals for: Test, Mr Patient

Referral Selection

Past Referrals | New Referral

Pathology | Radiology | Referral Letter

Addressee

☐ Free Text Addressee

Address

5. Search for the doctor, select their name and click **OK**.

Address Book Search

Name: Johnson ☐ Include Unavailable

Suburb:

Specialty: All

Type of Contact: All

Full Name	Work Phone Number	Address	Specialty	Notes	Work Fax Number	Mobile Phone	Home Phone	Email
Johnson, Dr Bob								

6. Once the correct doctor to send to is selected, You can type an appropriate subject for the letter. In this instance we have used **GP Referral Medical Objects Test Practice**.

Referrals for: Test, Mr Patient

Referral Selection

Past Referrals | New Referral

Pathology | Radiology | Referral Letter

Addressee

☐ Free Text Addressee

Copy To

Address

Letter Subject: GP Referral Medical Objects Test Practice

Other

Ref Date: 07/01/2016

Expiry Date: 07/01/2016

☐ No Renew Expected

☐ Send as PDF of letter

Today's Referrals

Address	Details	Problem
Dr Bob Johnson	GP Referral Medical Objects Test Practice	Undiagnosed

Current Encounter

Undiagnosed Problem

☐ Use Coded RFE

Referral: GP Referral Medical Objects Test Practice - Dr Bob Johnson

Once you have the letter subject, you can click the **Write Letter** button highlighted above.

7. From this screen you can type the letter.

Letter Editor

Testing Medical Centre
60 Park Street, Brisbane QLD 4000
Tel : 07 1234 5678

Dr Example Provider
M2355
000000000

7 January 2016

Dr Bob Johnson

Re:

Dear Dr Bob Johnson,

Re: Mr Patient Test DOB: 16/09/1967

Thank you for seeing Mr Patient Test, age 49 yrs, for opinion and management.

Presenting Problems:


Once you are ready, Click the highlighted **Distribute** button.




Zedmed will then ask you do you want to save this referral letter. Click **Yes**.

8. The Print/Send Screen will now appear. From here you need to make sure distribute now is selected. Also make sure that Send Electronically to all electronic recipients is ticked and then you can select the sending Doctor.

You can also print a copy of the document if necessary as well.

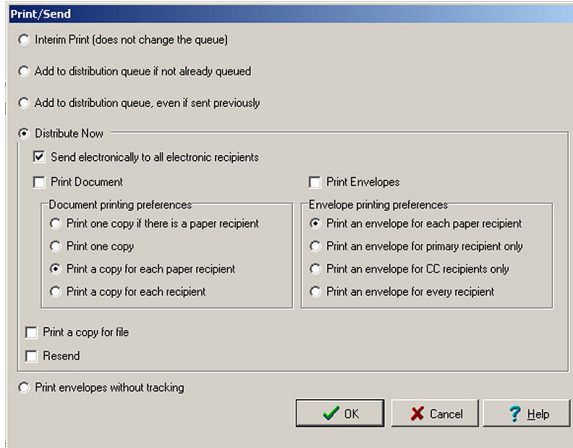

9. You can then click the **OK** button to send the report. Zedmed will then popup saying the following.

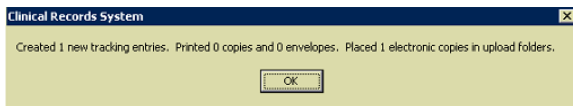


Sending Batch Letters via ZedMed

1. Click 'Batch Writer'



2. The list shown under 'Search Results' in the patients seen on these days

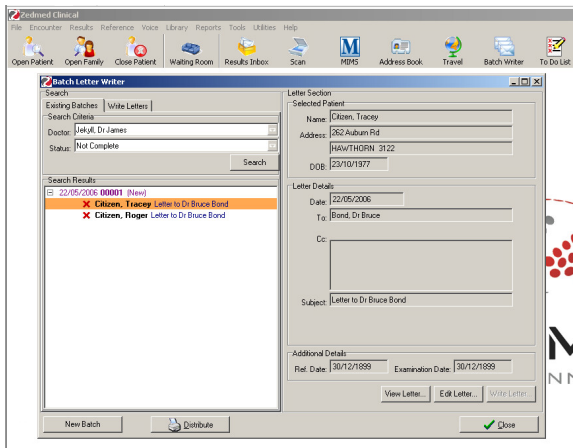

3. Change the tab to Write Letter to lookup a patient, click on the binoculars icon to search for a patient.

Batch Letter Writer

Search

Existing Batches: **Write Letters**

Search Criteria

Attendance Date From: 01/11/2012 To: 20/11/2017

6:00 AM 11:59 PM

Find Search

Selected Batch Details

Doctor: **Jekyll, Dr James** Batch No.: **00002**

Search Result

4. Type in patient's Surname in the 'Family Name' field. Now click 'find' to lookup the patient.

Find Patient

1. Family Name: Patient

Given Name

File #

DOB Gender

Medicare Number

Address

Home

Hospital UR No

Unique Patient Id

Unique Referral Id

☐ Search Archived and Deceased

2. **Find** New List (F4) Cancel Help

5. Now select the correct patient name, and then click 'Select'.

CRS Patient Selector

All Records for this Search Loaded

FILE #	Name	DOB	Address	PHONE	STATUS
12	Patient, Mr Test	12/12/1990			CUR

Other file numbers: New Include Archived and Deceased **Select** Cancel

5.

6. Now the patient has been selected, click the 'Write Letter' button to bring up the ZedMed letter writer.

Batch Letter Writer

Search

Existing Batches: Write Letters

Search Criteria

Attendance Date From: 01/11/2012 To: 20/11/2017

6:00 AM 11:59 PM

Find Search

Selected Batch Details

Doctor: **Jekyll, Dr James** Batch No.: **00002**

Search Result

20/11/2017 10:45am Patient, Mr Test

20/11/2017 10:43am Patient, Mr Test

20/11/2017 10:40am Patient, Mr Test

Letter Section

Selected Patient

Name: Patient, Mr Test

Address:

DOB: 12/12/1990

Letter Details

Date: 20/11/2017

To:

Co:

Subject: Letter

Template: Ref Letter - Standard

Summary:


Additional Details

Ref. Date: 20/11/2017 Examination Date: 20/11/2017

View Letter **Write Letter** Close

6.

7. From this stage you can follow the sending process as per the guide.

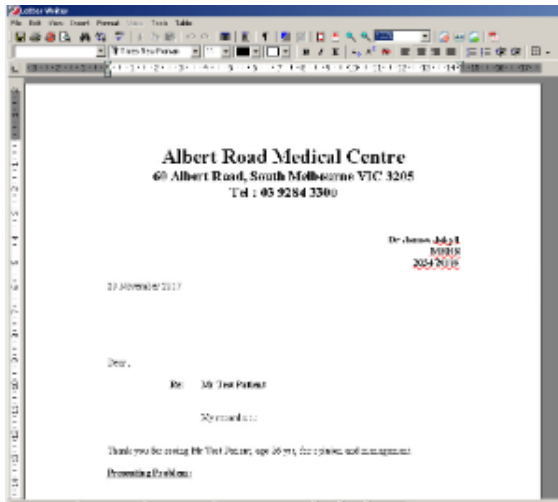

this will save time manually opening every patient if they usually type and send their results all at once.

Note:

 You can see a little envelope next to some names – This envelope tells the site what 'to doctors' (next to the patient's name) are setup to send electronically so they can easily send the ones that can go through.