


# Checking Result Delivery - Advanced

## Understanding the Delivery Report


The Explorer delivery report is capable of running many different types of reports. The default option is for checking on reports that have gone out for that particular day.

- 1. To check on the current day's reports, click **Reports**, and choose **View Document Delivery**.

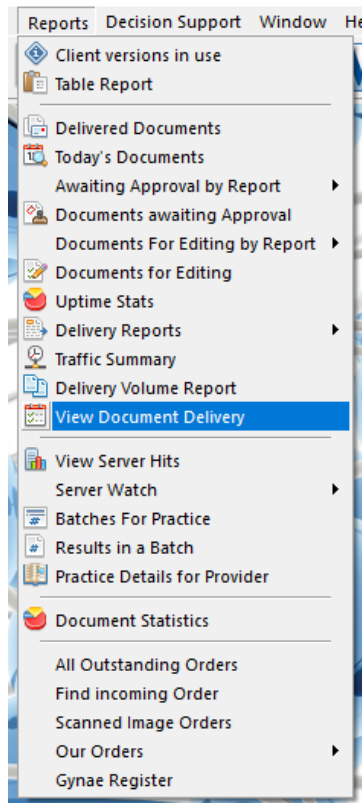

- 1 Understanding the Delivery Report
- 2 Reports that have been delivered
- 3 Locating a delivered report in referrer software.
- 4 Reports that have been delivered but not imported or consumed
- 5 Reports that have not been delivered
- 6 Troubleshooting Undelivered Reports
  - 6.1 Check Provider Number is registered
  - 6.2 Check Recipient site is currently online
- 7 Advanced Report Searching
  - 7.1 Search by Date Range
  - 7.2 Search by Patient Name
  - 7.3 Search by Doctor
  - 7.4 Combination Search
- 8 Re-Sending a Report

- 2. You will then be presented with a window like the one below. Simply click **OK** to view today's outgoing reports.

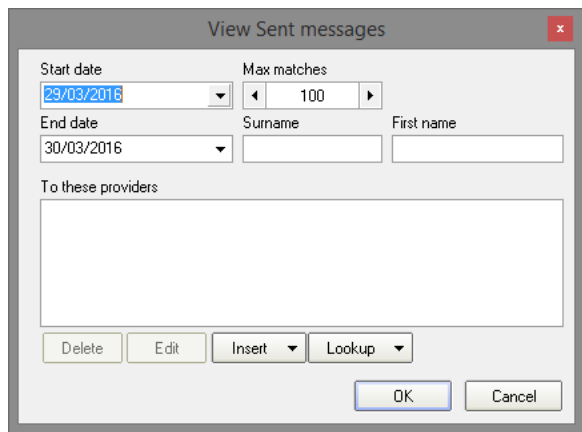

- 3. You will then be shown a list of reports that have been sent today.

Report Title	Delivered	Report Link/Copy For	Patient
CT Coronary Arteries	15.3.2016 08:05	<a href="#">Retrieve</a>	
MRI C Spine Cervical Radiculopathy (non Con)	15.3.2016 09:03	<a href="#">Retrieve</a>	
MRI C Spine Cervical Radiculopathy (non Con)	15.3.2016 08:08	<a href="#">Retrieve</a>	

The report is comprised of a number of fields. You may need to scroll to the right to see all fields depending on your screen resolution.

**Report Title** – e.g. X-ray Chest

**Delivered** – If populated, this shows the day and time that the report was downloaded by the recipient

**Report Link** – clicking this link will show you a copy of the report that was sent

**Copy For** – The name and provider number of the recipient

**Patient** – Patient name

**Observation Date** – Date and time that the test was performed.

**Delivery Type** – Shows the medical-Objects software type in use at the recipient site. They will be one of 3 types

Capricorn – PC Download Client  
CapricornCM – Mac Download Client  
EquatorDxTray – Medical-Objects clinical software

**Reviewed** – Only relevant when delivery type is EquatorDXTray

**Consumed** – Only populated when application acknowledgement is not created. Shows date and time that the report was picked up from the download location.

**Result Status** – Identifies if result was sent as a correction, interim result, or final result. Should show as Final in most cases.

**Application Ack** – There may be more than one application ack. One from our Capricorn to say the message has downloaded, and one from the receiving clinical system to say the report has been imported.

## Reports that have been delivered

When checking the report status, you first want to make sure the report has been delivered. So this would involve looking at the **Delivered** column and ensuring that there is a date and time in that column. If the entry is blank, then follow the troubleshooting steps listed [Here](#).

For this example, we will look at a number of files that have been delivered, and follow up on whether they have been imported at the recipient end.

In the previous picture there are 3 reports, all of which have a delivered date and time in the column. Let's look at the other fields to the right to see their status. See image below.

	Consumed	Result status	Application ACK	Application ACK
1	15.3.2016 08:13	Final result	CAPRICORN @ 15.3.2016 08:05	
2		Final result	CAPRICORN @ 15.3.2016 09:03	Genie v8.72 @ 15.3.2016 09:54
3		Final result	CAPRICORN @ 15.3.2016 08:08	MD3.15.3 @ 15.3.2016 08:16

You will see most reports have a 'Capricorn' entry in **Application Ack** column. This details the time that the download client (Capricorn) saved the report to a directory at the receiving end. If the report has not been delivered, there will be NO application ack entries.

As we know all 3 reports above have been delivered, we now look to see if they have been imported by the recipient's practice software.

You will either get an entry under the **Consumed** column, or under the second **Application Ack** Column. The differences between the two are explained as follows.

**Application Ack** – This is a receipt given to us by the recipient practice software, for example Medical Director, or Genie. It is essentially their software saying 'Thank you, file received' and details the date and time that it was imported by the clinical system. If this entry is populated, there is no doubt that the file has been imported by their software.

**Consumed** – **Consumed** means that the downloaded file is no longer in the location we delivered to. Indicating that the recipient's clinical software has picked the file up from that folder at that date and time. This is a more generic acknowledgement, created by Medical-Objects when the recipient system is not capable of giving us an **Application Ack**.

Let's look at the picture above, and explain each row by number.

1. The **Consumed** column has a date and time in there. This means that the report has been picked up from the download directory by the clinical system. So we now know that in this case, the file was downloaded by the Capricorn download client on 15/3 at 8:05am, and was picked up from the download location at 8.13am the same day.

2. The **consumed** column is blank, so we will move along and look in the 2nd **Application Ack** Column. It says in there '**Genie v8.72 @ 15.3.2016 09:54**'. In this case, we now know that the file was downloaded by the MAC Capricorn download client (Capricorn CM) on 15/3 at 9:03am, and then subsequently imported by Genie clinical system at 09:54am the same day.
3. The **consumed** column is again, blank, so we will move along and look in the 2nd **Application Ack** Column. It says in there '**MD 3.15.3 @ 15.3.2016 08:16**'. In this case, we now know that the file was downloaded by the Capricorn download client on 15/3 at 8:08am, and then subsequently imported by Medical Director clinical system at 08:16am the same day.

The only exception to above will be if you see the Delivery Type is displaying **EquatorDxTray**. This means they are using software similar to the Medical-Objects explorer as their clinical package. If you see **EquatorDxTray** as the delivery type, you should instead look under the **Reviewed** column for a date and time. This specifies exactly when they marked the report as reviewed. If this is blank, it means they have not yet reviewed the report.

See below for an example.

Delivery Type	Reviewed	Consumed	Result status	Application ACK	Application ACK
EQUATORDXTRAY	8.3.2016 13:57	N/A	Final result	EQUATORDXTRAY @ 8.3.2016 13:48	

## Locating a delivered report in referrer software.

If you receive a call from a referrer who claims they can't see a report, despite it appearing to be delivered and imported, the most common reason is that they are either looking in the wrong place in their clinical system, or they have inadvertently deleted the file. We suggest asking the referrer to;

- A) Check that the report is not in their holding file or doctor inbox
- B) Check that the report has not already been matched to the patient record.

In most cases when we follow up a 'missing' report, they are more often than not looking in the holding file when the report has already been matched to the patient record, or vice versa. The majority of these issues are often resolved by directing them to the correct spot.

If the site still cannot see the report, please notify (or have them call) the Medical-Objects helpdesk, and we will follow up with them.

## Reports that have been delivered but not imported or consumed

In some cases you may see something like the following.

Consumed	Result status	Application ACK	Application ACK
	Final result	CAPRICORN @ 15.3.2016 09:15	

There is a Capricorn ack to say the file has been downloaded, however there is no **consumed** ack, and no second **application** ack.

This indicates the file has been downloaded, but has not been imported by the clinical software on the receiving end. This could be for a number of reasons, such as

- a) The recipient clinical software could be on a scheduler to pick up the reports every hour or so. Try checking again after an hour has passed from the file downloading to see if a consumed or application ack is generated.
- b) The site prefers to manually import results at their own discretion and would have to perform an import for an acknowledgement to then be generated.
- c) The site could have changed clinical software and be looking in a different path for reports, in which case we need to fix the import paths.

If a site complains of not being able to view any incoming reports, despite them being downloaded at their end, one of the above conditions may be true. You can notify the Medical Objects helpdesk and we will follow up with the recipient site to resolve.

## Reports that have not been delivered

In some cases you may see that a report has no entry in the delivered column (see example below). This means the report has not been downloaded by the recipient's download client

US Abdominal ☒ **Retrieve** JOHNSON, BOB (005562AW)

There can be a few reasons why a site has not downloaded reports. The most common 2 reasons are:

1. The sites download client is offline or having problems. This can be because of Internet being down at the site, download client not running, server or computer replacement, disk corruption, or similar reasons.
2. The provider number that you are sending to is not registered on the Medical-Objects network.

Below are some steps to troubleshooting undelivered reports.

## Troubleshooting Undelivered Reports

If you see a report does not have a delivered entry in the column, you can trouble shoot by checking the following in order.

1. Is the provider number you have sent to registered with Medical-Objects?
2. Is the recipient site currently online?

Below are instructions for checking each of these conditions.

Check Provider Number is registered

Medical-Objects deliver reports using the recipient provider number as our addressing method. It is crucial to report delivery that the reports are addressed to a provider number that is registered on our system. The most common reason for undelivered reports is that they are addressed to an unregistered provider number: therefore our system does not know where to deliver them to.

To check if the provider number you are sending to is online with Medical-Objects, open your explorer client and choose **Lookup Provider** from the **Providers** Menu

Providers View Navigate Tools  
Lookup Provider

Enter the provider number you are sending to. If you get 'no matches' then this means the number is not in our system.

☒ Australian Directory
 ☐ Local Providers
 ☐ Online Providers Only

Surname <space> Firstname <Return>
 Max 500
 Search Type
 Provider Number


005565TW

Name	Practice Name	Phone
No matches		

It is advised to check again by searching the provider by surname, as we may have the provider connected at the required location but with a different provider number. Set the search type to 'Name' and search by surname <Space> firstname as below.

☒ Australian Directory
 ☐ Local Providers
 ☐ Online Providers Only
 ☐ Show Inactive Providers



Sunname <space> Firstname <Return>
 Max 500
 Search Type Name
Add...

Name	Practice Name	Phone	Address	Provider ID
 JOHNSON, BOB	Testing Practice	(07) 12345678	TESTER 123 Test Street	005562AW

In the example above, if the provider shows up and the practice is correct, you should change the provider number you have in your in house system to the one that is registered with us and then re-send the message from your in house system (RIS/LIS/Other)

If you get no results, or the provider is not registered at the location you want to send to them at, we would require the referrer to email us a request to add the provider to their download client. Once this is done, any messages already in the system should then get delivered. You should not have to resend to that provider number.

Beware of entries in the directory that do not have the routable icon next to their names. If the routable icon is not there, then this means that the provider is in our directory, but NOT connected to any download client. This can happen when someone is added to the system for an installation which has not occurred yet. See below for example.

Routable Icon		<table><thead><tr><th>Name</th><th>Practice Name</th><th>Phone</th><th>Address</th></tr></thead><tbody><tr><td>CAMPBELL, SAM</td><td>MO Testing Capricorn</td><td></td><td>MAROOCHYDORE</td></tr><tr><td>CAMPBELL, SAM</td><td>MOA Test Capricorn</td><td>(07) 54566000</td><td>MAROOCHYDORE 102 Villes Road</td></tr></tbody></table>	Name	Practice Name	Phone	Address	CAMPBELL, SAM	MO Testing Capricorn		MAROOCHYDORE	CAMPBELL, SAM	MOA Test Capricorn	(07) 54566000	MAROOCHYDORE 102 Villes Road
Name	Practice Name	Phone	Address											
CAMPBELL, SAM	MO Testing Capricorn		MAROOCHYDORE											
CAMPBELL, SAM	MOA Test Capricorn	(07) 54566000	MAROOCHYDORE 102 Villes Road											
Unroutable														

## Check Recipient site is currently online

If a recipient's download client is not online, then your reports will queue in the system until such time that their download client is available again.

If you would like to view the current status of the recipient's download client, open your [Delivery Report](#) and click the the [Provider Lookup](#) section. Then follow the [Provider Lookup](#) instructions on the [Delivery Report](#) page to run the search. Next, find the provider number of the recipient that you are wanting to send to and make sure that the provider number is currently valid by looking for the green tick in the online column.

Name	Practice	Phone	Address	Provider ID	Online
MO Testing Capricorn	MO Testing Capricorn	(07) 54566000	MAROOCHYDORE 102 Villes Road	4555	✓
MOA Test Capricorn	MOA Test Capricorn	(07) 54566000	MAROOCHYDORE 102 Villes Road	4556	✗

If there is red cross in the online column, please inform the Medical-Object helpdesk.

If you have confirmed the provider number is connected, and the site is online, yet reports are still not going through, then this would indicate another issue, most commonly that the file path at the site where we download reports to has been changed or deleted. You should also notify the helpdesk of this to enable us to rectify the issue.

## Advanced Report Searching

Sometimes you will receive an enquiry from a practice, who is chasing up a single report, or who perhaps haven't received any reports within a particular time frame, such as the past 2 weeks.

Rather than search through a massive list of sent reports trying to find the one you are looking for, you can narrow down your search results to return only the ones you are interested in.

### Search by Date Range

When you run the delivery report, you are presented with the View Sent Messages search box. The default setting (clicking OK) will return reports sent for the current day only. However, if you wish to see all reports sent within a set period, you can alter the Start and End dates to show all reports sent between those dates. In the example below this search will show all reports sent for the month of February.

View Sent messages

Start date	Max matches	
1/02/2016	100	
End date	Surname	First name
1/03/2016		

### Search by Patient Name

If a practice is looking for a report for a particular patient, the easiest way to find this report is to set a date range when you know the report would have been sent between (as above), and then put in the

patient surname. The search example below will return all reports for patients with a surname 'Smith' that were sent in the month of February. If you want to further narrow the search, you can also put the patient's first name.

View Sent messages

Start date: 1/02/2016 Max matches: 100

End date: 1/03/2016 Surname: Smith First name:

## Search by Doctor

If a practice says that a particular doctor has not received reports, you can also search for all reports sent to that doctor within a date range. To do so, Enter your date range, and choose lookup, and then lookup staff records from the menu.

Start date: 1/02/2016 Max matches: 100

End date: 1/03/2016 Surname: First name:

To these providers:

Delete Edit Insert Lookup

Lookup Staff Records Ins

Search for the doctor you want to view reports for, and click on their name.

☐ Australian Directory ☒ Local Providers ☐ Online Providers Only ☐ Show Inactive Providers

Surname <space> Firstname <Return> johnson bob Max 500 Search Type Name Add...

Name	Practice Name	Phone	Address
JOHNSON, BOB	Testing Practice	(07) 12345678	TESTER, 123 Test Street

This will then allow you to see all reports sent to that particular doctor within that specified date range, and should look similar to below.

View Sent messages

Start date: 1/02/2016 Max matches: 100

End date: 1/03/2016 Surname: First name:

To these providers:

Dr Bob JOHNSON (005562AW)

Tip:



To see all reports sent to a particular practice, use the method above, but search by 'organisation name' and lookup the practice by name. You can then hold the CTRL key to mouseover and select multiple doctors.

## Combination Search

Now that you are familiar with the search methods above, you can combine all 3. The example below shows how you can search for all reports, sent in February, for patients with surname 'Smith', that were sent to Dr Bob Johnson.

**View Sent messages**

Start date: 1/02/2016 Max matches: 100

End date: 1/03/2016 Surname: Smith First name:

To these providers:

Dr Bob JOHNSON (005562AW)

Delete Edit Insert Lookup

## Re-Sending a Report

If you need to re send an already delivered report to another doctor, you can do this from within the Medical-Objects explorer.

First, find the report you wish to re-send using the search methods above.

Click the **retrieve** link on the delivery report.

Report Title	Delivered	Report Link	Copy Fo
Consultation Note (C POWELL)	17.3.2016 10:57	<a href="#">Retrieve</a>	LEUNG, MI

This will load the body of the report in the viewer window.

**Test PATIENT**

Address: 15/102 Wises Road MAROOCHYDORE QLD 4558 Phone: Not recorded

DOB: 02/04/1984 0478-4843-9239-FB0693102DA0 Report Date: 17/03/2016

Referral by: ASSOC PROF MICHAEL LEUNG CC

Retrieve All Results Print Result Create Response Provider Note

**Consultation Note (C POWELL)**

Re: Test Patient  
DOB: 12/12/12

15/102 Wises Road  
Maroochydore, QLD 4558  
Australia

Hi there,

This is a test message from Medical-Objects to see if your connection is working for delivery of results.

You may delete this message after 2 days.

Kind Regards,  
Cameron Powell  
Medical-Objects Pty Ltd  
11281 / www.medical-objects.com.au  
Web: http://web.medical-objects.com.au/  
Phone: (07) 5456 6000  
Fax: (07) 5421 0220

"Australia's Fastest Secure Messaging"

Now, click on **Documents**, from the menu, and choose **Send a Copy to Another Doctor**

**Documents** Create Document Provide

Format Print

Unreviewed Results

Documents for Editing

**Send a Copy to Another Doctor**

Show Document Audit Log

Search for the recipient by name or provider number, and click on their name

**Lookup Provider**

☐ Australian Directory ☒ Local Providers ☐ Online Providers Only ☐

Surname <space> Firstname <Return>

johns bob

Name	Job Title	Practice Name
JOHNSON, BOB		Testing Practice

The report will now have been sent. You should check your delivery report to ensure it has been delivered.

If you have any questions or any issues with the Explorer delivery report, Please contact the Medical-Objects helpdesk on (07) 5456 6000.